


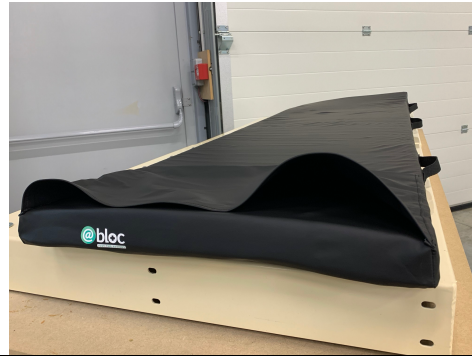

Instructions for use of the MEDEVAC BRAN002ARM @bloc stretcher mattress

Version v4 February 2025

The MEDEVAC stretcher mattress BRAN002ARM incorporates foams, resulting from tests that have allowed us to find the best physiological compatibility to maximize the comfort and safety of the patient, considerably reducing the risk of pressure ulcers. They are intended to be positioned on MEDEVAC® stretchers to optimize patient comfort.

Restrictions on Use

- The user must ensure that the mattress is suitable for the stretcher.
- This mattress is intended for users trained in patient positioning and the risk of pressure ulcers.
- Do not wash our cushions and mattresses in the washing and disinfection cabin.
- Never use abrasive cleaners or scraper sponges.
- Do not attach sticky strips to the mattress but if necessary to the sides of the stretcher, always making sure to insert a drape to avoid contact of the adhesive on the mattress or cushion.
- If an alcohol-based product for disinfecting the patient's skin is used, it is imperative for safety reasons (risk of inflammation and/or maceration), to ensure that it does not accumulate under the patient.
- Prolonged use of alcohol or nail polish remover may impair the conductivity and elasticity of the mattress or cushion covering and impair its waterproofness and/or cause it to shrink.
- Make sure the patient is not allergic to polyurethane.

Characteristics

- The mattress is made of high-density mouse (40kg/m³) and viscoelastic foam (85Kg/m³)
- Composition of the upper outer shell: Polyurethane (polyurethane-coated polyamide fabrics)
- Composition of the lower outer shell: PVC (PVC coated cotton backing)
- Traditional stitching assembly process
- The mattress is equipped with a pocket for a head cushion
- The mattress is equipped with 3 handles on each side allowing the passage of the support straps. These handles can also be used for transporting the mattress
- A strap is provided to keep the mattress in the rolled position during transport
- A waterproof carrying bag is available as an option (reference BRANSAC02)

Dimensions :

- Length 1900 mm, width 560 mm and thickness 50 mm


Mattress rolling

To facilitate its transport, the mattress can be rolled. To do this:

- Roll on the floor while compressing the mattress.
- Belt the mattress with the strap provided

For unwinding:

Proceed in reverse.

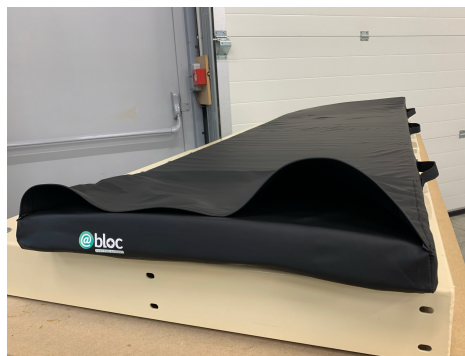

Washing, disinfection and maintenance

- Mattresses or cushions can be cleaned by wiping with a soft cloth dampened with water containing a neutral detergent or with a ready-to-use solution (e.g. ANIOS Surfa'safe, ALFAKARM Septalkan, HUCKERT'S Umonium 38 Medical Spray...) or soaked wipes (e.g. ANIOS Oxy Wipes; Wip'ANIOS Excel; HUCKERT'S Umonium 38 Medical Tissues Laboratory...).
- In general, surface disinfectants containing a set of active ingredients such as aldehyde, quaternary compounds, guanidine derivatives and active oxygen, have proven their effectiveness. But always choose those formulated specifically for equipment, less film-forming and aggressive.
- In situ cleaning and disinfection can also be carried out on the coating with hot water by hand and a neutral detergent or with a sodium hypochlorite solution (0.1% or 1000 parts per million of available chlorine).
- Products such as Ercesolvent©, Remove©, GILBERT non-stick,© etc. used to remove adhesive residues (Elastoplast©, sticky drapes), should be thoroughly rinsed and dried immediately after use. Failure to do so may result in the buildup of reagent damaging the polyurethane coating. All products based on acetone, ether and other volatile solvents must be used with the same precautions.
- It is essential that items dry thoroughly after cleaning. Thus, we should never cover our mattresses with drapes, gels, heated mattress toppers, or put them back on the stretcher before they are completely dry.
- Finally, it is imperative for hygiene reasons to replace the mattress or cushion as soon as it is damaged (tear, cut, cracks, burns, etc.).

Special precautions

Before any use, it is necessary to ensure that the mattress is in good condition and that it is properly installed.
The outer coating of the mattress is not intended for permanent contact with the patient's skin.
The mattress can be rolled up and strapped in for easy transport. For prolonged storage, it is preferable to put the mattress horizontally.
The use of this mattress is not a substitute for the monitoring necessary to prevent pressure ulcers. Pressure ulcer prevention requires identifying risk factors and planning position changes every 2 to 3 hours, or even more frequently.

Special conditions of storage and use

Min temperature	Storage : -15 °	Use : + 15°
Max temperature	Storage : + 45 °	Use : +25°
Hygrometry (%HR) min	Storage : 30%	Use : 30%
Hygrometry (%HR) max	Storage : 70%	Use : 70%


End of life

At the end of its life, the mattress must be cleaned and disinfected normally, then placed in the system intended for the disposal of waste of the general type.



Tel : + 33(0) 2 41 17 49 49
Fax: + 33(0) 2 72 22 11 43
info@abloc.eu

8 rue des Frères Montgolfier
49240 Avrillé
France

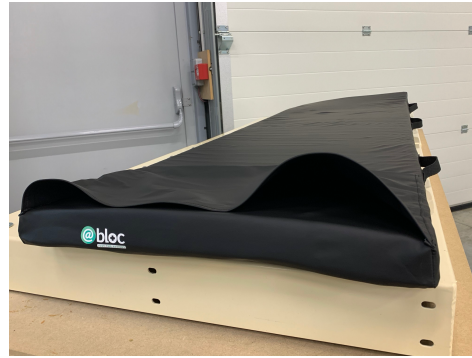

Normative references

Manufacturer: @bloc
Complies with European regulations
CE Class 1 Medical Device
Inflammability
BS7175:1989 crib5 for top fabric
M1 for the lower fabric
Tissue in contact with the patient: Antistatic according to ISO 2878

Any serious incident that occurs in connection with one of our mattresses must be notified to @bloc and the ANSM in accordance with the regulations in force.

For any further information, please contact us.

